



Slim WISC face shield assembly

Components

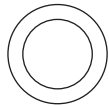

roll of double-sided tape




clear plastic sheet for shield visor



Elastic strip




6 inch x 1 inch foam

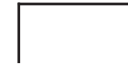


staples and stapler


Assembly

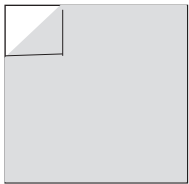

1. Peel film off both sides of plastic if needed – not all face shield visors have plastic film.

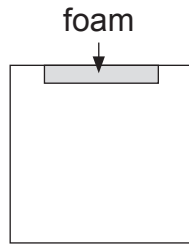

2. Put double-sided tape on foam and center foam along top edge of shield visor. Push foam down firmly so it sticks.

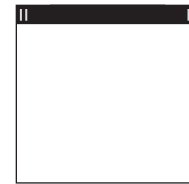

3. Double staple elastic on both sides close to edge to finish shield.

(elastic should cover foam, but do not staple through foam)

This is a quick reference guide. Point your phone camera at this QR code to watch a video of assembly.

